



HOKKAIDO UNIVERSITY

Title	An instrument-free, screen-printed paper microfluidic device that enables bio and chemical sensing
Author(s)	Mohammadi, Saeed; Maeki, Masatoshi; Mohamadi, Reza M. et al.
Citation	Analyst, 140(19), 6493-6499 https://doi.org/10.1039/c5an00909j
Issue Date	2015
Doc URL	https://hdl.handle.net/2115/60219
Rights(URL)	https://creativecommons.org/licenses/by/3.0/
Type	journal article
File Information	Supplementary information.pdf



Supporting Information

Instrument-Free, Screen-Printed Paper Microfluidic Device That Enables Bio and Chemical Sensing

Saeed Mohammadi ^a, Masatoshi Maeki ^b, Reza M. Mohamadi ^c, Akihiko Ishida ^b, Hirofumi Tani ^b,
Manabu Tokeshi ^{b,d,e,f} [†]

^a Graduate School of Chemical Sciences and Engineering, Hokkaido University, Kita 13 Nishi8, Kita-ku, Sapporo 060-8628, Japan

^b Division of Biotechnology and Macromolecular Chemistry, Faculty of Engineering, Hokkaido University, Kita 13 Nishi 8, Kita-ku, Sapporo 060-8628, Japan

^c Department of Pharmaceutical Science, Leslie Dan Faculty of Pharmacy, University of Toronto, Toronto, Ontario M5S 3M2, Canada

^d ImPACT Research Centre for Innovative Nanobiodevices, Nagoya University, Furo-cho Chikusa-ku, Nagoya 464-8603, Japan

^e Innovative Research Centre for Preventive Medical Engineering, Nagoya University, Furo-cho Chikusa-ku, Nagoya 464-8603, Japan

^f Institute of Innovation for Future Society, Nagoya University, Furo-cho, Chikusa-ku, Nagoya 464-8603, Japan

[†] *Corresponding author: TEL: +81-11-706-6744, E-mail: tokeshi@eng.hokudai.ac.jp*

Optimisation of pipetting volume

Because assay reagents are expensive, optimization of pipetting volume of solution into paper sample zones is important. Moreover, using a large volume of solution than the sample zones can contain can lead to cross contamination as a secondary problem. On the other hand, using a smaller volume of solution than the sample zones can contain can lead to miss-pipetting in multistep assays. 0.01 M fluorescein solution was traced using a fluorescence microscope as a model to investigate the optimum pipetting volume of solution in the sample zones. The investigation procedures included two simple steps: first different volumes of fluorescein solution (0.3-0.7 μL) were spotted in the sample zones and second, the fluorescence images of the sample zones were obtained with the fluorescence microscope. As shown in Figure S-1, 0.2, 0.3, and 0.4 μL volumes of the solution were not sufficient to fill all the sample zones. In addition for the reasons mentioned above, 0.6 and 0.7 μL volumes and more were not desirable.

From these results, we considered 0.5 μL as the optimum pipetting volume in a sample zone (Figure S-1).

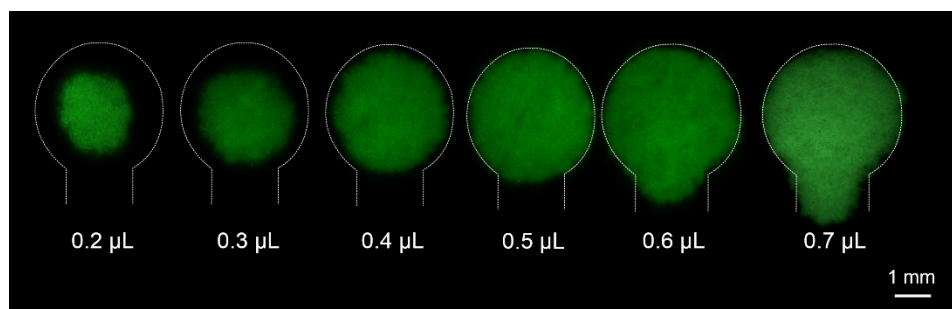


Figure S-1. Optimisation of pipetting volume in sample zones and sufficiency of solution immobilization.